

The Synd-EZ kit provides a solution for damage to the syndesmosis joint which connects the tibia and fibula bones of the lower leg. This simple yet effective device is designed to help stabilize the ankle and can be used with or without a fracture plate.

Features & Benefits

Versatile design

- Ideal for use with or without a fracture plate
- Washer included for use without a fracture plate
- Top hat is low profile and self centering on fracture plates

Adjustable knotless construct


- Simplifies surgical procedure

Counter traction suture & tension handle

- Integrated for predictably smooth tension of construct

Synd-EZ					
Part #	Width of Button (mm)	Length of Button (mm)	Description	Length of Loop (mm)	QTY/Box
11223	3.2	12	Synd-EZ Titanium kit	-	1

Synd-EZ Instrumentation					
Part #	Diameter (mm)	Joint	Description	Sterile/Unsterile	Single-use/Reusable
11272	-	Small Joint	Synd-EZ drill guide	Unsterile	Reusable

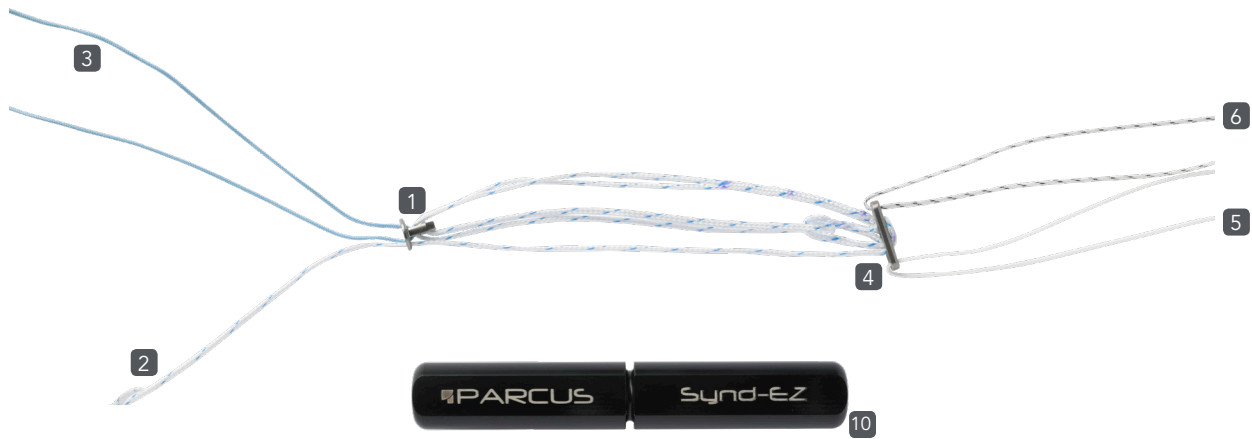

Synd-EZ Kit Components	
1	Top hat
2	White/blue adjustment suture
3	Solid blue counter-traction suture
4	Button
5	Solid white passing suture
6	White/black toggle suture
7	3.5mm cannulated drill bit
8	2.4mm passing pin
9	1.5mm K-wire
10	Black tension handle
11	Washer

